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Digital Automático de Pulso

N ML
266
INMETRO 2009

Z-43

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Características:

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Z-43

Instruction Manual

- Features
- Automatic
- Dual Person Memory Zones
- 50 Memories for each zone
- Average of last 3 Readings
- Date/time indications
- Easy Operation

**CONSULT YOUR PHYSICIAN IF THERE IS ANY VARIATION IN
YOUR READINGS. THE USE OF THIS PRODUCT DOES NOT
DISPENSE THE MEDICAL ADVICE.**

DO NOT USE THE PRODUCT ON INVASIVE MEASURING

INTRODUCTION

Congratulations on your purchase of the new blood pressure monitor wrist type model **Z-43**. This fully automatic, wrist mounted blood pressure monitor is both easy to use and ideally suited for daily measurements. The enhanced and comprehensive features will make **Z-43** the only blood pressure monitor you will need.

With this blood pressure monitor, you can store up to **50 Measurements per Memory Zone**, and it is ideal for users who want to monitor and track the blood pressure on a regular basis. This wrist watch monitor is compact and portable, making it ideal for home use and travel.

Attention: Please read this manual carefully before use. For specific information on your own blood pressure, contact your physician. Please be sure to keep this manual.

ABOUT BLOOD PRESSURE

What is Blood Pressure?

Blood pressure is the pressure exerted on the artery tube while blood flows through the arteries. The pressure measured when the heart contracts and sends blood out of the heart is systolic (highest) blood pressure. The pressure measured when the heart dilates with blood flowing back into the heart is called diastolic (lowest) blood pressure.

Why Measure Your Blood Pressure?


Among the various health problems afflicting modern people, problems associated with high blood pressure are by far the most common. High blood pressure dangerously correlates with cardiovascular diseases and high morbidity. Therefore blood pressure monitoring is a necessity to identify these risks.

BLOOD PRESSURE STANDARD

The World Health Organization (WHO) has developed a blood pressure standard, according to which areas of low and high risk blood pressure are identified. This standard, however, is a general guideline as individual's blood pressure varies among different people and different age groups, etc.

It is important that you consult with your physician regularly. Your physician will tell you your normal blood pressure range as well as the point at which you will be considered at risk.


For reliable monitoring and reference of blood pressure, keeping long-term records is recommended.

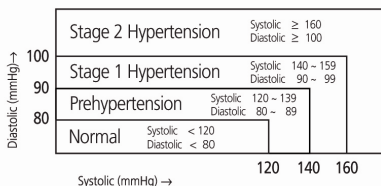BLOOD PRESSURE FLUCTUATION

Blood pressure fluctuates all the time!

You should not be over worried if you encountered two or three measurements at high levels. Blood pressure changes over the month and even throughout the day. It is also influenced by season and temperature.

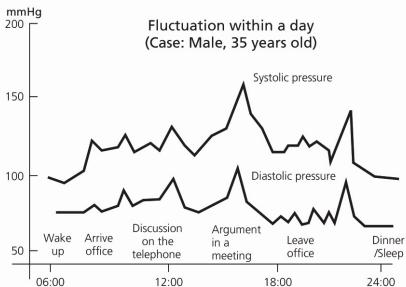
HOW Z-43 WORKS


Z-43 uses oscillometric method to detect your blood pressure. Before the cuff starts inflating, the device will establish a baseline cuff pressure, which is equivalent to the air pressure. After the cuff inflates to block your blood in the artery, the deflation process starts. During the deflation of the wrist cuff, the monitor is detecting the pressure oscillations generated by the beat-to-beat pulsatile. Any muscle movement during this period of time will cause measurement error. After detecting the amplitude and the slope of the pressure oscillations during the deflation process, your **Z-43** will determine the systolic and diastolic pressures for you, and your pulse rate is detected at the same time.

CARRYING CASE

Please store the blood pressure monitor in the carrying case as shown below:

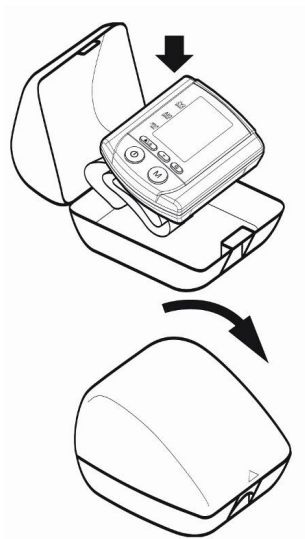NAME/FUNCTION OF EACH PART

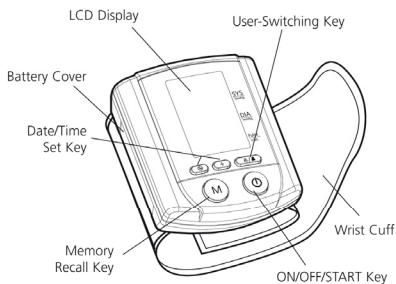
Accessories for unit :

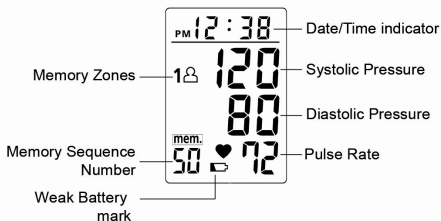2 "AAA" size, 1.5V alkaline batteries included in the carrying case with the unit.

DISPLAY EXPLANATIONS

Display:

ICONS:

	Memory zone 1: Appears when the monitor is operated in memory zone 1.
	Memory Zone 2: Appears when the monitor is operated in memory zone 2.
	Weak Battery Mark: Appears when batteries should be replaced.
	Pulse Mark: Shows the pulse rate per minute.
mem.	Memory Mark: Shows the number of stored measurements.
	Measurement Error: Exceeding measurement range. Measure again quietly. Wrap the cuff correctly and keep wrist steady during measurement. If error keeps occurring, contact the customer service center.
	Air Circuit Abnormality: Measure again quietly. If error keeps occurring, contact the customer service center.
	Pressure Exceeding 300 mmHg: Switch the unit off and measure again quietly. If error keeps occurring, contact the customer service center.
	Data Error: Remove and reload the batteries. If error keeps occurring, contact the customer service center.


INSTALLING BATTERIES

1. Battery Type: 2 Alkaline LR03 (AAA) size

2. Place your thumb on the side of the battery cover. Pull the battery cover in the direction of the cuff to open.

3. Install the batteries following (+) and (-) marks.

4. Replace the cover by placing the end tabs of the battery cover in the slots on the top of the battery compartment. Click in the other end to secure the battery cover.

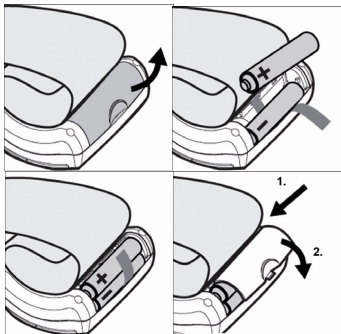Replace the batteries if:

1. The weak battery mark appears in the display.

2. Nothing appears in the display when the power is switched on.

As the supplied batteries are for test only, they may be discharged earlier than batteries you buy in stores. Replace the batteries in pairs. Remove batteries when unit is not in use for extended periods of time.


Batteries are hazardous waste. Do not dispose them together with the household garbage.

APPLYING THE CUFF

1. Remove the watch, jewelry, etc. prior to attaching the wrist monitor. Clothing sleeves should be rolled up and the cuff should be wrapped on bare skin for correct measurements.

2. Apply cuff to left wrist with palm facing up as Fig. A.

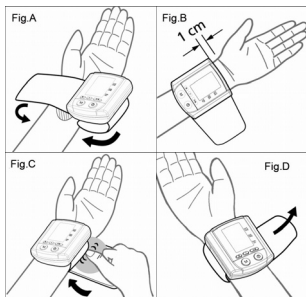3. Make sure the edge of the cuff is about 1 cm from the palm as Fig. B.

4. In order to ensure accurate measurements, **fasten the Velcro strap securely** around your wrist so there is no extra space between the cuff and the wrist as Fig. C.

If the cuff is not wrapped tight enough, the measurement values may be incorrect.

5. If your physician has diagnosed you with poor circulation on your left arm, carefully place the cuff around your right wrist as shown in Fig. D.

CORRECT MEASURING POSTURE

1. Place your elbow on a table so that the cuff is at the same level as your heart as Fig. A


Note:

Your heart is located slightly below your armpit a bit to the left of the middle of your chest. Relax your entire body, especially the area between your elbow and fingers.

2. If the cuff is not at the same level as your heart or if you cannot keep your arm completely still throughout the reading, use a soft object such as a **folded towel** to support your arm as Fig. B. Do not allow hard objects to come in contact with the wrist cuff.

3. Turn your palm upwards.

4. Sit upright on a chair, and take 5-6 deep breaths. Avoid leaning back while the measurement is being taken as Fig.C.

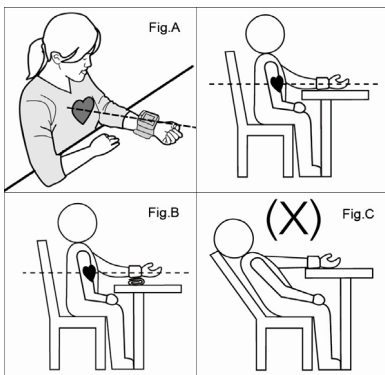
HELPFUL TIPS

Here are a few helpful tips to help you to obtain more accurate readings:


- Do not measure your blood pressure immediately after consuming a large meal. To obtain more accurate readings, please wait one hour before measuring.
- Do not smoke or drink alcohol before measuring your blood pressure.
- You should not be physically tired or exhausted while taking a measurement.
- It is important that you relax during measurement. Try to take a 15 minute rest before taking measurements.
- Do not take measurements if you are under stress or under tension.
- Take your blood pressure at normal body temperature. If you are feeling cold or hot, wait a while before taking a measurement.
- If the monitor is stored at very low temperature (near freezing), have it placed at a warm location for at least one hour before using it.
- Wait about 5 minutes before taking the next pressure measurement.

MEASUREMENT PROCEDURES


1. Press the **User-Switching key** to select Memory Zone 1 or Memory Zone 2.

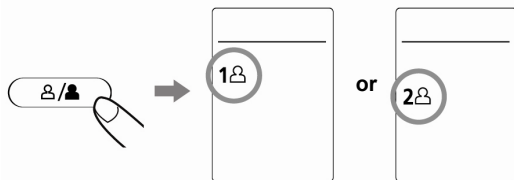
2. Place the cuff on the wrist. Press the  **key**. All digits will light up, checking the display functions. The checking procedure will be completed after about 2 seconds.

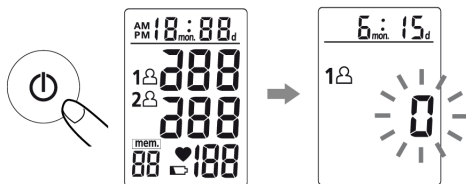3. After all symbols appear, the display will show a blinking “0”. The monitor is **“Ready to Measure”** and will **automatically** inflate the cuff to start measurement.

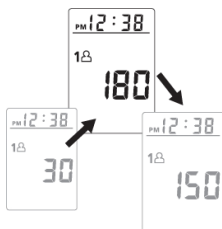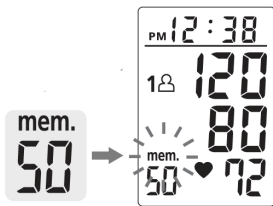4. When the measurement is completed, systolic pressure, diastolic pressure and pulse will be shown simultaneously on the LCD screen. **The measurement is** then automatically stored into the pre-designated memory zone.

5. This monitor will re-inflate automatically to **approximately 220 mmHg** if the system detects that your body needs more pressure to measure your blood pressure.

Note:


1. This monitor automatically switches off approximately **1 minute** after last key operation.
2. To interrupt the measurement, simply press the  key. The cuff will deflate immediately after  key is pressed.
3. During the measurement, do not talk or move your arm or hand.

RECALLING VALUES FROM MEMORY

1. The Monitor has two Memory Zones (1 and 2). **Each zone can store up to 50 measurements.**

2. To read memory values from a selected Memory Zone, use the **“User-Switching Key”** to select a Memory Zone (1 or 2) from which you want to recall values from. Then, press the **“M”** key, and the values will be shown on the LCD display with an assigned memory sequence number.

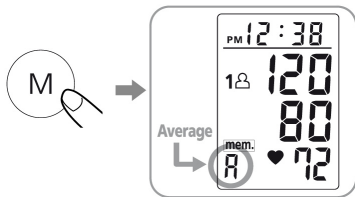3. The monitor will display an average of the previous measurements first. Every new press of the **MEMORY** key calls for one set of prior values stored. Latest memory will be recalled first.

The calculated average will be applied to a maximum of the last 3 memories.

CLEARING VALUES FROM MEMORY


1. Press the **User-Switching** key to select the memory zone 1 or memory zone 2.

2. Continue to press the “M” key for approximately 5 seconds, then the data in the pre-designated memory zone can be erased automatically.

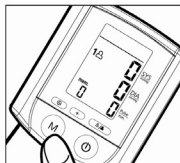**Press and hold
for 5 seconds.**


Note:

The data in both memory zones (zone 1 and zone 2) can be erased if any of the batteries is moved away. (Suggest the users to first record the last readings before replace or move away the batteries.)

TIME ADJUSTMENT

1. To adjust the date/time, press . The display will show a blinking number showing the month, for example.
2. Change the month by pressing the  **key**. Each press will increase the number by one in a cycling manner. Press  again to confirm the entry, and the screen will show a blinking number representing the date.
3. Change the date, the hour, and the minute as described in Step 2 above, using the  **Key** to change and  to confirm the entries.
4. A "0" will reappear as the Blood Pressure Monitor is ready for measurement again.

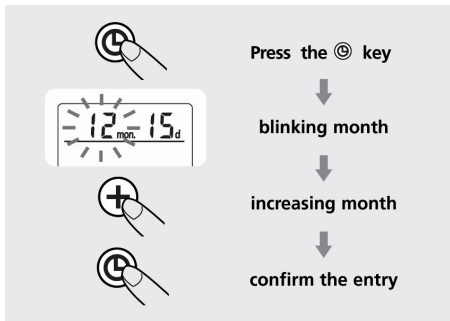TROUBLESHOOTING

If any abnormality will arise during use, please check the following points.

Symptoms	Check Points	Correction
No display when the  key is pressed	Have the batteries run down?	Replace them with two new alkaline batteries.
	Have the batteries polarities been positioned incorrectly?	Re-insert the batteries in the correct positions.
EE mark shown on display or the blood pressure value is displayed excessively low or high.	Is the cuff placed correctly?	Wrap the cuff properly so that it is positioned correctly.
	Did you talk or move during measurement?	Measure again quietly. Keep wrist steady during measurement.
	Did you shake the wrist with the cuff on?	

Note: If the unit still does not work, contact the customer service center. Under no circumstance should you disassemble and repair the unit by yourself.

CAUTIONARY NOTES

- 1.** The unit contains high-precision assemblies. Therefore, avoid extreme temperatures, humidity, and direct sunlight. Avoid dropping or strongly shocking the main unit, and protect it from dust.
- 2.** Clean the blood pressure monitor body and the cuff carefully with a slightly damp, soft cloth. Do not press. Do not bend the pre-formed cuff inside out. Do not wash the cuff or use chemical cleaner on it. Never use thinner, alcohol or petrol (gasoline) as cleaner.
- 3.** Leaky batteries can damage the unit. Remove the batteries when the unit is not used for a long time.
- 4.** Do not press the ON/OFF/START key when the cuff is not placed around the wrist.
- 5.** The unit should not be operated by children to avoid hazardous situations.
- 6.** Do not disassemble the main unit or cuff, nor use any kind of tool to open it.
- 7.** If the unit is stored near freezing point, allow it to acclimate at room temperature before use.
- 8.** For users diagnosed with common arrhythmia (atrial or ventricular premature beats or atrial fibrillation), diabetes, poor circulation of blood, kidney problems, or for users suffered from stroke, or for unconscious users, the device may have difficulty in determining the proper blood pressure.
- 9.** To stop operation at any time, press the ON/OFF/START key, and the air in the cuff will be rapidly exhausted.
- 10.** Once the inflation reaches 300 mmHg, the unit will start deflating rapidly for safety reason.
- 11.** It is advisable to calibrate your blood pressure monitor once a year.

SPECIFICATIONS

Power Source:	3V DC, Two LR03 (AAA) Alkaline Batteries
Measurement method:	Oscillometric
Measurement range :	Pressure : 40~250mmHg; Pulse: 40~199 beats/minute
Pressure sensor :	Semi-conductor
Accuracy :	Pressure: ± 3 mmHg; Pulse: $\pm 5\%$ of reading
Inflation:	Pump driven
Deflation :	Automatic Pressure release valve
Memory capacity:	50 mem. for each zone x2 zones
Auto-shut-off :	1 minute after last key operation
Operation enviroment:	10°C~40°C (50°F~104°F); 40%~85% RH max.
Storage enviroment:	-10°C~60°C (14°F~140°F); 10%~95% RH max.
Dimensions :	81 (L) X 69 (W) X 30 (H) mm
Weight :	120 g (without batteries)
Wrist circumference:	Adulto: 13.5~22 cm (5,3"~8,7")
Limited users:	Adult users

Type BF

Device and cuff are designed to provide special protection against electrical shocks.

*Specifications are subject to change without notice.

WARRANTY

1. The main unit is guaranteed for 2 years and the cuff, being a fragile part, its guarantee is stipulated for 3 months from date of acquisition and this guarantee is effective only by presenting the **RECEIPT OF PURCHASE** stating the buyer's name, date of purchase, name, and address of the distributor.
2. Only defective products with manufacturing problems are covered by this guarantee. Devices with malfunctioning problems caused by misuse or improper handling even if not intentionally, are not covered by guarantee.
3. This guarantee does not cover damages by dropping down the device, direct exposure to high humidity, water, battery leakage, and the batteries included in the device, repairs executed by unauthorized person, violation of the device by the customer or by unauthorized person.
4. The guarantee does not cover the cost of round-trip transport for repair, the result or damage caused by improper use of the device and other expenses not specified here.
5. The guarantee is good only in countries where this product is officially commercialized by the manufacturer or the local distributor and authorized by the manufacturer, always with the presentation of the receipt of the purchase.
6. This is a non-transferable guarantee; it means it is good for the first consumer only.
7. During the period of this warranty the manufacturer will repair the product after 30 days from the date of receiving the device for repair.
8. The manufacturer reserves the right to replace the defective product by a new one in case of need and this is decided by the manufacturer.
9. All repairs executed during the period of warranty do not extend the deadline of the guarantee.

10. After the guarantee is expired, the manufacturer has to manage to reserve components and to supply technical assistance during 5 years which will be paid by the customers.

11. During the warranty period all services of maintenance and of consumed components under bad manipulation are not acceptable.

12. The periodical revision/calibration of the product is not covered by the guarantee.

13. The technical assistance and warranty are restricted to the Brazilian territory only.

14. This guarantee does not cover damages by misuse or improper handling, even if not intentional, dropping down the device, direct exposure to high humidity, water, and sunshine.

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Other regions: 0800-555-750.

Other countries: please contact your local distributor.

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